
Caremaxx®

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Infrarood lamp • Lampe infrarouge
Infrared lamp • Infrarot Lampe

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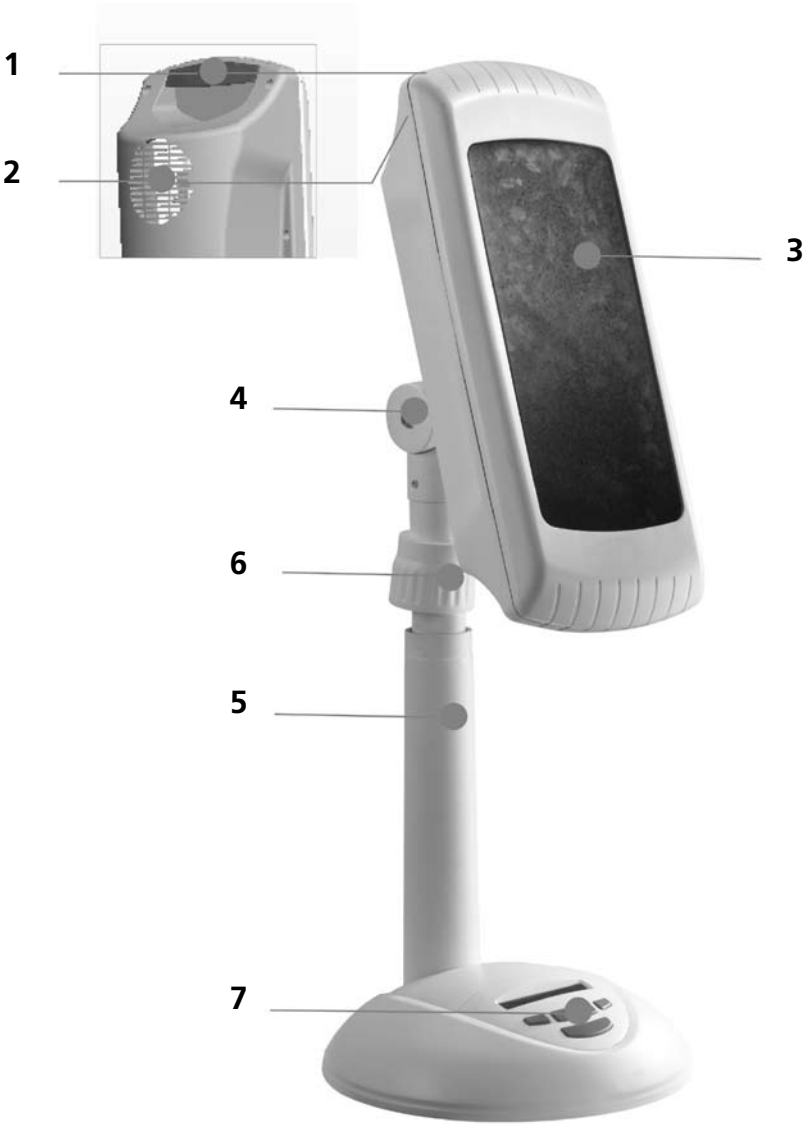1. Handvat
2. Luchttoevoer
3. Lichtfilter
4. Scharnier
5. Voet
6. Verstelring
7. Bedieningspaneel

Description générale

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E-mail : info@caremaxx.nl
Site Internet : www.caremaxx.nl

Introduction

Thank you for purchasing a Caremaxx infrared lamp.

This infrared lamp is a Caremaxx quality product.

Please make sure to carefully read this user manual before using the instrument. Keep the user manual in a safe place in case you wish to consult it later on.

General description

1. Handle
2. Air supply
3. Light filter
4. Hinge
5. Foot
6. Adjusting ring
7. Operating panel

Safety precautions

- Use the product only according to the instructions in the user manual.
- Inappropriate use will invalidate the warranty.
- The product is not meant for commercial use.
- Before connecting the instrument, make sure to check whether the voltage mentioned on the instrument corresponds to the local mains voltage.
- Repair jobs may only be performed by authorized service workshops.
- Infrared radiation may heat up any objects in the vicinity of the lamp so they become quite hot.
- Be careful – this poses a potential hazard for burns! Keep the instrument out of reach of babies, children, and mentally handicapped persons.
- Watch the instrument well, never leave it alone as long as it is switched on.
- Never touch the instrument, and especially the light filter ☺, if the instrument is running, as the lamp will heat up! Also, do not touch it for the first 15 minutes after switching off the instrument.
- Using the instrument in the vicinity of flammable, explosive or heat-sensitive objects or gases is strictly prohibited. Keep the instrument away from flammable objects (such as textiles).
- Never use the instrument in a dusty environment.
- Never cover up either the instrument or the vent holes. Keep at least 35 cm free space around the instrument so as to prevent overheating.
- Keep the instrument away from other heat sources and avoid direct exposure to the sun.
- Never use the instrument at room temperatures below 10°C or above 35°C.
- Do not expose the instrument to strong impact.
- Never attach the instrument to a wall or a ceiling.
- Pull the plug out of the socket before cleaning the lamp.
- Allow the lamp to cool off for at least 15 minutes before cleaning it or putting it away.
- Keep this instrument away from water.
- Never immerse any part of the instrument in water or any other liquid.

- Do not use the instrument in humid rooms or out in the open.
- You may clean the instrument with a soft, slightly moist cloth soaked with a small amount of a neutral detergent, then dry the products thoroughly. Do not plug in the lamp unless and until it is completely dry.
- Pull the plug out of the socket if the instrument is not in use. Keep the instrument – preferably in its original packaging - in a dry and dust-free place.
- The lamp's infrared radiation is unsuitable for tanning.
- Do not stare straight into the lamp; instead, keep your eyes closed to prevent damaging them.
- Never use the instrument if the light filter ③ is damaged, broken, or missing.
- Protect the treated body part from chilling.
- Do not fall asleep during treatments.
- Never use the instrument if you are on pain killers—which are prone to affect your heat sensitivity. Consult your physician before using the instrument.
- Never use the instrument if you have any swellings or inflammations.

Delivery and packaging

The packaging contains:

- 1 infrared lamp (650Watt)
- 1 user manual

The packaging material may be either reused or recycled. Make sure that used packaging materials are sent to the pertinent waste management facility.

Any transport damages detected while unpacking should be reported your supplier immediately.

Application

The Caremaxx infrared lamp may be used at home to treat muscle problems and to alleviate pain. Infrared radiation penetrates deep into the body, which is therefore heated from within. Heating the cutaneous tissue will widen the blood vessels and enhance the blood flow and the metabolism, thus also accelerating the evacuation of waste products from the body.

The Caremaxx infrared lamp lends itself, among others, for the treatment of:

- stiff muscles
- muscle pain
- common colds
- (lower) back pain
- rheumatic problems

There should be some noticeable improvement after 7 sessions; otherwise please consult your family physician.

Always consult a medical doctor before using the instrument if you are running a fever, while you are using pain killers, if you are suffering from (infrared) light hypersensitivity, swellings, inflammations, diabetes, heart problems, multiple sclerosis, thromboses, or coagulation disorders.

Operation

Place the lamp on a stable, level surface (for instance a table). Set the desired height by

using the adjusting ring ⑥. Firmly tighten the adjusting ring to prevent the lamp from descending during treatment. While adjusting the height, hold the lamp by the hinge ④ or the handle ①.

Completely unwind the power cable. Insert the plug into the socket. Turn the instrument on by pressing the main switch (ON/OFF). The green LED lamp will light up. The lamp is turned off by once more pressing the main switch (ON/OFF).

The desired time may be set by pressing the "Time +" and "Time -" buttons (between 0 and 30 minutes). The desired light intensity (high, medium, low) is set by adjusting the mode button. The medium position is the pre-programmed standard setting that is automatically activated when switching on the lamp. This time period may be adjusted during treatment. The lamp automatically switches off as soon as this time period expires. Disconnect the plug from the socket and allow the lamp to cool off before putting it away.

The treatment may be interrupted by using the time-out button. The lamp will switch off, but the remaining time period stays in the display, which makes it possible to proceed with the treatment a little later by pressing the time-out button once more.

The instrument should never be used more than 15 minutes at a time. It is preferable to perform several short treatments instead of one longer treatment.

The tilt angle of the lamp may be adjusted

as well. Hold the lamp by its handle ⑤ and position it as desired, making sure not to touch the light filter and burning yourself.

Keep at least 50-60 cm distance from the lamp to prevent the skin from overheating. If the heat feels uncomfortable to you, simply increase the distance until the heat feels comfortable.

Cleaning and maintenance

Always be sure to disconnect the plug from the socket and allow the instrument to cool off before cleaning it. The outside may be cleaned with a moist cloth. Make sure that no water can enter the instrument. Never use any abrasive cleaners, abrasive sponges, or liquids such as alcohol, benzene, or acetone. After cleaning, dry the instrument thoroughly with a dry cloth.

Defective halogen lamps are to be replaced by an authorized workshop.

The halogen lamp is not covered by the warranty. You are not entitled to free halogen lamp replacement.

Waste management

Neither this instrument nor the lamps are to be disposed of with the household waste. Consumers are obliged to send all electric or electronic instruments, whether or not they contain any hazardous materials, to a depository for environmental pollutants in their home town or to their specialized dealership to ensure their environmentally sound disposal.

If you have any questions regarding

waste management, please contact your municipality or your specialized vendor.

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