

Caremaxx®

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Lite Pad lichttherapie • Lite Pad lumineothérapie
Lite Pad light therapy • Lite Pad Lichttherapie

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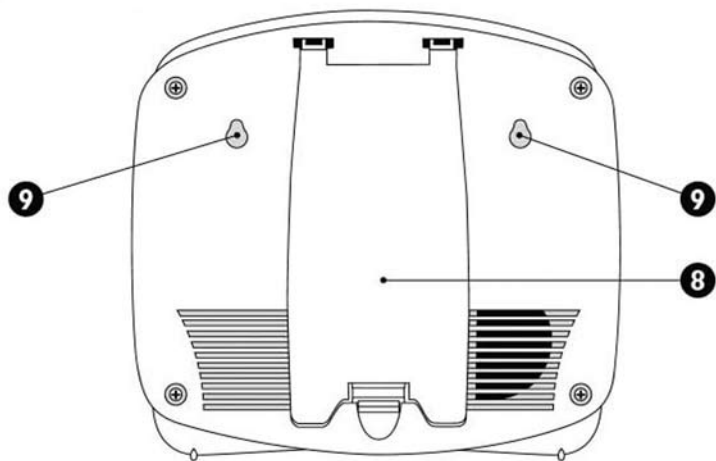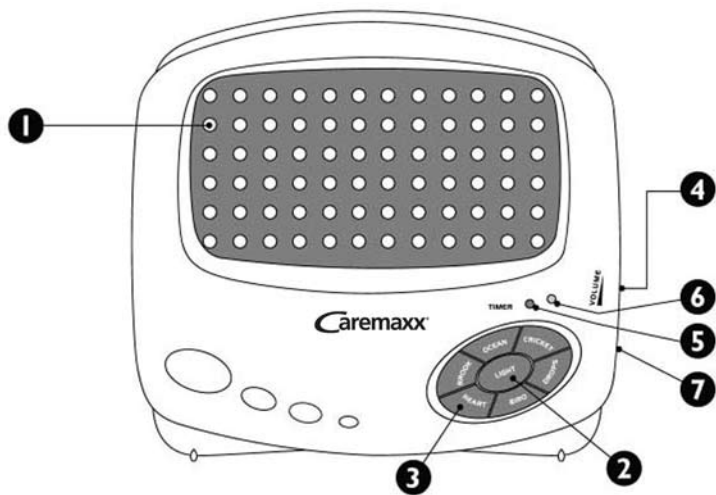E-mail : info@caremaxx.nl
Site Internet : www.caremaxx.nl

1 Introduction

Thank you for purchasing a Caremaxx Lite Pad.

This Lite Pad is a **Caremaxx** quality product.

Some of us have problems during the short days of fall and winter with their lack of adequate day light. This may be compensated with the Lite Pad.

Make sure to carefully read this user manual before using the instrument. Keep the user manual in a safe place in case you wish to consult it later on.

2 Important

- Use the Lite Pad only according to the instructions in the user manual.
- Inappropriate use will invalidate the warranty.
- The Lite Pad is not meant for commercial use.
- Before connecting the instrument, make sure that the voltage mentioned on the instrument corresponds to the local mains voltage.
- Repair jobs may only be performed by authorized service workshops.
- Keep this instrument away from water; never use this instrument in humid rooms.
- Never immerse the Lite Pad in water or any other liquid – neither partially nor completely.
- Keep the instrument out of reach of children.
- This instrument ought to be used in regular lighting conditions and not in a completely dark room, due to the difference in contrast.

- If you are suffering from depression or symptoms of winter depression, please consult your physician.
- The Lite Pad is not a medical device.
- The Lite Pad cannot be used as a replacement for medical therapy.

3 Delivery and packaging

The packaging contains:

- 1 Lite Pad
- 1 Adapter
- 1 user manual

The packaging material may be either reused or recycled. Make sure that the used packaging material is sent to the pertinent waste management facility in charge. Any transport damages detected while unpacking should be reported your supplier immediately.

4 Winter depression

Winter depression (or seasonal affective disorder) is sometimes also referred to as winter blues. This disorder occurs mostly in the wintertime.

Primary symptoms include:

Sleep disorders:

Extreme fatigue, sleep deprivation, difficulties falling asleep at night and phase misalignment of the sleep/wake cycle.

Eating disorders:

In the wintertime, patients with winter depression tend to eat extra-high amounts of carbohydrates (such as bread or pasta) as well as sugar (for instance in the form of chocolate), which may cause them to put on quite some weight. Weight gains of 10 kilograms

are not exceptional. This additional weight is then lost again in the spring.

Depressive, dull or listless mood without perceptible cause.

Patients with winter depression are much more prone to lose their balance in the wintertime, they may react irritably to their environment and ultimately isolate themselves. They often lose interest in social interaction and may suffer from concentration problems. These symptoms of depression will vanish again once the days are getting longer.

5 Light therapy

The emphasis in treating winter depression is on the amount of light the patient is exposed to each day. This is known as light therapy. In this context the patient spends a certain amount of time daily in front of a strong light source. The lamp should not be too far removed from the eyes. The best time of day is in the morning because the lamp makes it easier to wake up properly.

The light intensity of the **Caremaxx** Lite Pad reaches a maximum of 5,000 lux and is thus many times higher than that of a regular lamp. In the wintertime, outdoor light intensities vary between 2,000 and 10,000 lux. Indoor light intensities may be even lower and vary between 0 and 500 lux. A sunny summer day, on the other hand, easily reaches a light intensity of 100,000 lux. Thus, regular lamps are inadequate, as they are normally below 500 lux.

Using the Lite Pad for 15 to 30 minutes on a daily basis is sufficient to reduce

complaints related to winter depression. Do not stare directly into the light. Lite Pad light therapy is a simple and effective way to compensate for the shortage of light without using medication.

6 How to operate the instrument

The duration of each session depends on the distance between the face and the light source. The smaller the distance between the face and the light source, the briefer the session.

You may use the Lite Pad while doing something else, such as reading, writing, telephoning, etc. Keep in mind, however, that you will have to remain close to the lamp during treatment.

Place the Lite Pad on a solid surface (e.g. table) using the foldable stand.

Put the plug into the mains socket and switch on the device using the On/Off switch (2).

Assume position between 30 and 60 cm from the lamp.

We advise you not to look into the light.

Once the light is on, you can activate the timer by pressing 1x on the On/Off switch (2). The timer will flash 1x to show that the timer has been configured to 15 minutes.

If you press the On/Off button 2x, the timer will be configured to 30 minutes. The timer will flash 2x to show that the timer has been configured to 30 minutes.

You can turn off the device by pressing down the On/Off switch (2) for a few seconds.

Follow a regime of 5 to 7 consecutive days.

There is no need to complete the entire treatment in a single session. If it suits you better, you may easily interrupt the treatment and resume it later that day. After treatment, turn off the instrument by pressing the On/Off switch (2) and disconnect the plug from the socket.

The first few times you are using the Lite Pad you may feel some eye fatigue or a slight headache. These side effects usually vanish quite soon.

Make sure to consult your physician before using the Lite Pad:

- If you are or have been suffering from severe depression;
- If you are hypersensitive to light (for instance because of diabetes or epilepsy);
- If you are suffering from an eye disease;
- If you are using certain prescription drugs (especially antidepressants, psychotropic drugs, or malaria tablets);
- In case of high blood pressure;
- If you have any questions regarding light therapy in view of your specific health status.

7 Natural sounds

You are able to choose from 6 relaxing natural sounds (3):

1. Ocean
2. Crickets
3. Drops
4. Bird
5. Heart
6. Brook

Press the appropriate button to play the sound of your choice.

You can activate the timer by pressing the same button again. 15 minutes are added to the timer every time you press the button:

Press 1x > timer 15 minutes, timer indicator (5) flashes 1x

Press 2x > timer 30 minutes, timer indicator (5) flashes 2x

Press 3x > timer 45 minutes, timer indicator (5) flashes 3x

Press 4x > timer 60 minutes, timer indicator (5) flashes 4x

You can set the volume of the sound by using the volume button (4).

8 Cleaning and maintenance

Always be sure to disconnect the plug from the socket and allow the instrument to cool down before cleaning it. The outside may be cleaned with a moist cloth. Make sure that no water can enter the instrument. Never use any abrasive cleaners, abrasive sponges, or liquids such as alcohol, benzene, or acetone.

9 Waste management

Neither this instrument nor the lamps are to be disposed of with the household waste. Consumers are obliged to send all electric or electronic instruments, whether or not they contain any hazardous materials, to a depository for environmental pollutants

in their home town or to their specialized dealership to ensure their environmentally sound disposal.

If you have any questions regarding waste management, please contact your municipality or your specialized vendor.

10 Warranty and conditions for repair

In cases where you want to use your warranty, you must always contact your dealer or our customer service. If you need to send the device to our customer service, you must include the complaint and a copy of your receipt.

The following conditions apply to the warranty:

1. For **Caremaxx** products, a warranty period applies of three years after the purchase date. This can be proven by means of the receipt or invoice.
2. All complaints resulting from defective material or manufacturing faults, will be repaired for free within the warranty period.
3. Cases where warranty applies, do not automatically lead to extension of the warranty period, neither for the device itself, nor for the parts that have been replaced.
4. Excluded from the warranty are:
 - a. All damages caused by incompetent handling, e.g. not correctly following the user manual.
 - b. Damages caused by repairs performed by the buyer or another unauthorized person.
 - c. Transport damages, caused during

transport from the seller to the user or when sending the device to the customer service.

- d. Accessories that are subject to wear.
5. The manufacturer does not accept any liability for direct or indirect consequential damage caused by the device. This also applies, if the damage to the device is recognized under this warranty.

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